


# DISPOSABLE TROCAR

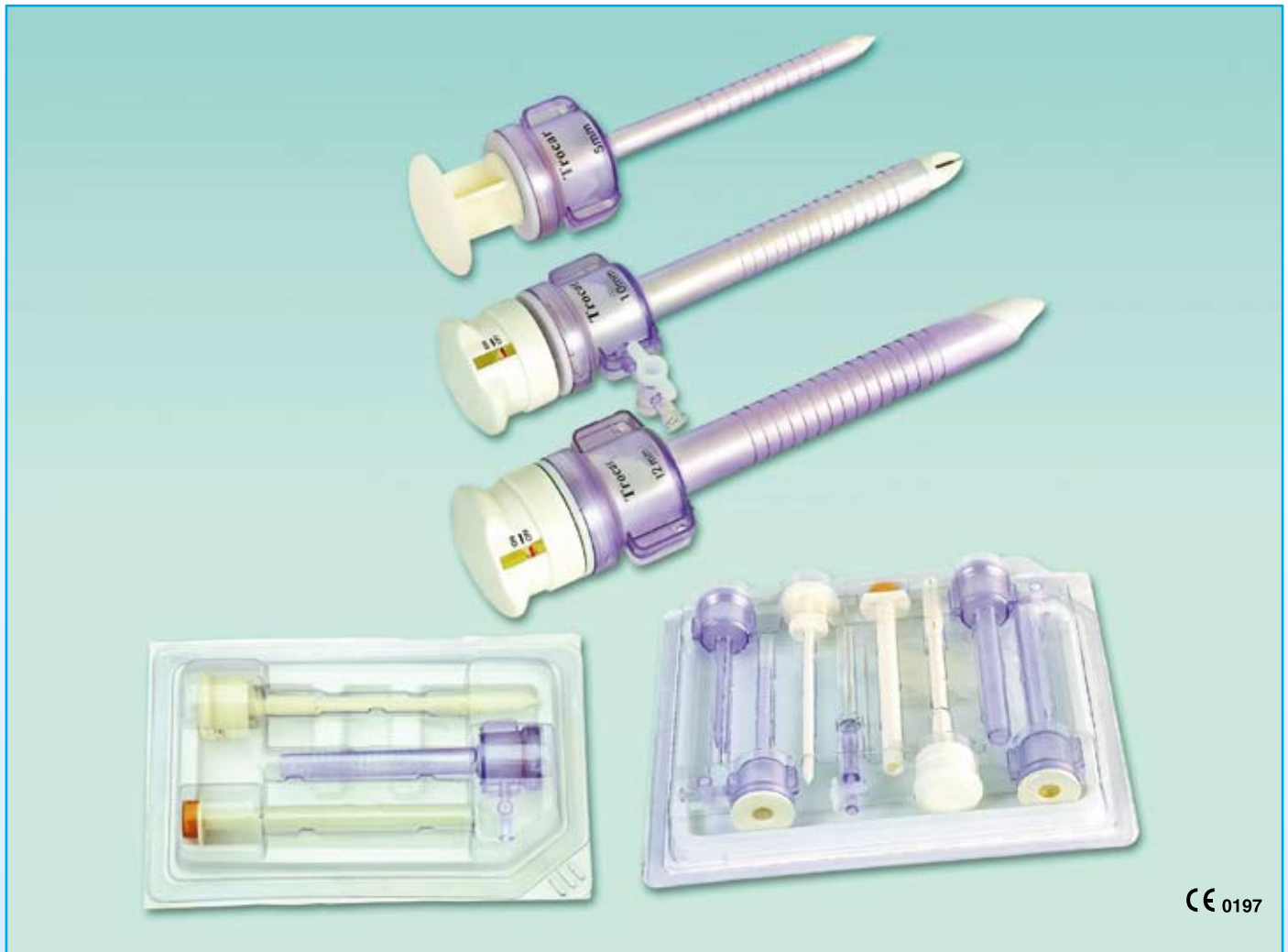

## The trocar is an important instrument in endoscopic surgery.

Three sizes are available: **5 / 10 / 12 mm**.

The Stapleline trocar set consists of 2 sizes of trocar sleeves and obturators, together with adapter sleeves, Veress needle and endobag. Endobags are available in two sizes: medium and large.

Its properties are as follows:

1. The Stapleline trocar is a sterile single-use instrument.
2. The 10 mm trocar has a spring-mounted protective cover, which protects the scalpel.
3. The trocar is marked, so that the surgeon can judge the position of perforation protection. During the puncture process, the blade of the scalpel only protrudes 2.5 mm: this prevents any injuries to inner organs.

4. The Stapleline trocar sleeve is gas-tight and has a CO2 insufflation port.
5. The Stapleline trocar sleeve has grooves which hold the trocar in the tissue, at the same time they stabilize the sleeve which thus has a smaller diameter. In this way, the wound can be kept as small as possible. There are 2 lateral openings to fasten the trocar sleeve during surgery.

### Ordering Information

| Trocar Set 1             | Order number: |
|--------------------------|---------------|
| 2 x Trocar sleeve 5 mm   | C04-02-Z05B   |
| 2 x Trocar sleeve 10 mm  |               |
| 1 x Trocar insert 5 mm   |               |
| 1 x Trocar insert 10 mm  |               |
| 1 x Reducer              |               |
| 1 x Veress needle 120 mm |               |
| 1 x Endobag 10 x 13,5 cm |               |

| Trocar Set 2             | Order number: |
|--------------------------|---------------|
| 2 x Trocar sleeve 5 mm   | C04-02-Z06B   |
| 1 x Trocar sleeve 10 mm  |               |
| 1 x Trocar sleeve 12 mm  |               |
| 1 x Trocar insert 5 mm   |               |
| 1 x Trocar insert 10 mm  |               |
| 1 x Trocar insert 12 mm  |               |
| 1 x Reducer              |               |
| 1 x Veress needle 120 mm |               |
| 1 x Endobag 10 x 13,5 cm |               |

| Separate Items             | Order number:  |
|----------------------------|----------------|
| 10 x Reducer               | B04-10-1004    |
| 10 x Veress needles 120 mm | B04-10-1005    |
| 10 x Endobags 12 x 15 cm   | B04-10-1006-01 |
| 10 x Endobags 15 x 22 cm   | B04-10-1006-03 |



CE 0197

## Der Trokar ist ein wichtiges Instrument in der endoskopischen Chirurgie.

Es sind drei Größen erhältlich: **5 / 10 / 12 mm**.

Das Stapleline Trokar Set besteht aus zwei Größen von Trokarhülsen und Obturatoren, sowie aus Reduzierhülsen, Veress Needle und Endobag. Endobags können in zwei Größen ausgewählt werden: mittel und groß.

Dieses neue Trokar Set ist sehr sicher und operationspraktikabel. Seine Eigenschaften sind wie folgt:

1. Der Stapleline Trokar ist ein steriles Einweginstrument.
2. Der 10 mm Trokar besitzt einen federnd gelagerten Obturator, welcher das Skalpell schützt. Der Trokar ist mit einer Markierung versehen, die es dem Arzt beim Operieren ermöglicht, die Position des Perforationsschutzes zu beurteilen. Beim Punktieren kommt die Schneide des Skalpells nur 2,5 mm hervor. Damit wird eine Verletzung der inneren Organe vermieden.

3. Die Stapleline Trokarhülse ist gasdicht und verfügt über einen CO2 Insufflation Port.
4. Die Stapleline Trokarhülse hat Rillen, die den Trokar im Gewebe halten, gleichzeitig stabilisieren sie die Hülse, die dadurch einen kleineren Durchmesser erhält. Dadurch kann die Inzision möglichst klein gehalten werden.
5. Es gibt zwei seitliche Öffnungen zur Befestigung der Trokarhülse während der Operation.

### Bestell Informationen

| Trokar Set 1             | Bestellnummer |
|--------------------------|---------------|
| 2 x Trokarhülse 5 mm     | C04-02-Z05B   |
| 2 x Trokarhülse 10 mm    |               |
| 1 x Einführung 5 mm      |               |
| 1 x Einführung 10 mm     |               |
| 1 x Reduzierer           |               |
| 1 x Veress needle 120 mm |               |
| 1 x Endobag 10 x 13,5 cm |               |

| Trokar Set 2             | Bestellnummer |
|--------------------------|---------------|
| 2 x Trokarhülse 5 mm     | C04-02-Z06B   |
| 1 x Trokarhülse 10 mm    |               |
| 1 x Trokarhülse 12 mm    |               |
| 1 x Einführung 5 mm      |               |
| 1 x Einführung 10 mm     |               |
| 1 x Einführung 12 mm     |               |
| 1 x Reduzierer           |               |
| 1 x Veress needle 120 mm |               |
| 1 x Endobag 10 x 13,5 cm |               |

| Artikel separat            | Bestellnummer  |
|----------------------------|----------------|
| 10 x Reduzierer            | B04-10-1004    |
| 10 x Veress needles 120 mm | B04-10-1005    |
| 10 x Endobags 12 x 15 cm   | B04-10-1006-01 |
| 10 x Endobags 15 x 22 cm   | B04-10-1006-03 |